


FlexRx™ – 28-Cavity Weekly Standard Multi-Dose Medication Management System

Streamlined Efficiency for Pharmacy Teams

Next generation FlexRx™ – 28-Cavity Weekly Standard seamlessly merges time-saving innovation with sustainability. Tailored for efficient pharmacy teams, it boasts a unique design for quick label application and incorporates our sustainable Bio-PET blister.

Key Features & Benefits

- 1** Raised alignment posts on the blister ensure precise positioning for smooth and effortless label application
- 2** New medically approved Bio-PET blister contains 33% plant-based PET (available with and without hooks)
- 3** Blister corners absorb energy when medication is dropped during filling, decreasing the risk of "jumping" medications
- 4** 28 high-capacity blister cavities per package to use as a 1-week package with 4 medication passes/day
- 5** Hooks allow filled blisters to hang for convenient storage or in a medication cart
- 6** Standard A4 label (8.5" x 11.69") assures compatibility with all laser printers
- 7** Summary grid provides guidance during filling and allows for easy verification of filled product
- 8** Perforated label on all cavities allows easy access to medication
- 9** Printing on every cell indicates patient name, medication enclosed and the date/time it should be taken with shaded columns to visually cue the time of day to take medication
- 10** Each medication pass time is perforated and can be detached for "on-the-go" portability, creating an easy-to-use medication cup with push-and-tear access to medications

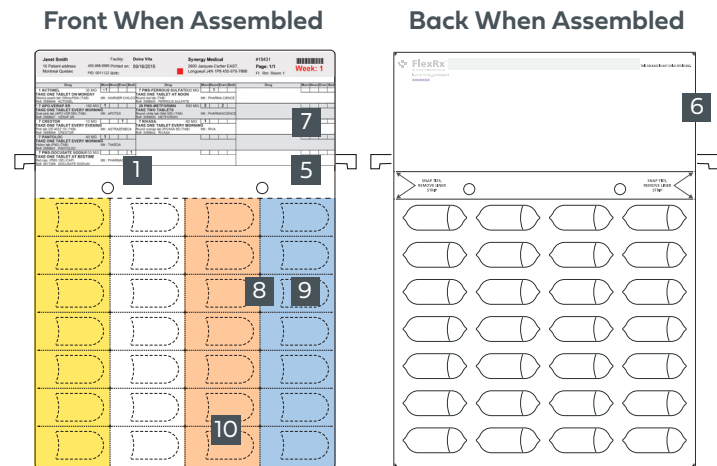

Compatible with BD Parata XF™, ULTRA™, Assist™, VBM 200F / TruCard™ and BlistMED / ATDPS DOB 208



TRY IT OUT WITH OUR STARTER KIT
TINTROKIT-FLEXRX2

DESCRIPTION	MCKESSON #	JONES ITEM #	CASE QUANTITY
FlexRx™ – 28-Cavity Weekly Short Label	193495	TLAS_24069-001	500
FlexRx™ – 28-Cavity Weekly Blister with Hooks – BioPET	193493	TBLIS_23968-00	500
FlexRx™ – 28-Cavity Weekly Blister without Hooks – BioPET	202026	TBLIS_24346-00	500
FlexRx™ – 28-Cavity Weekly Sealing Tray		TTRAY_24507-00	1
6" Hard-Seam Roller		T12210	1
Starter Kit – FlexRx™ – 28-Cavity Weekly		TINTROKIT-FLEXRX2	1

FlexRx™ – 28-Cavity Weekly Standard Multi-Dose Medication Management System

Sealing Instructions



1 Place blister: Begin by placing your empty blister on the sealing tray. Ensure accurate filling of the blister.



2 Prepare label: Identify the corresponding printed label and remove the perforated liner strip by peeling with the tab from the front.



3 Align label: After removing the perforated strip, align the label using the raised alignment posts as guides and press it firmly onto the blister.



4 Apply label: Fold the label to access the liner, peel off, and smooth the label to secure it to the blister.



5 Ensure secure adhesion: Use a roller to apply even pressure to the label in a vertical direction to ensure a secure seal.



6 Repeat process: Repeat the process for the remaining blisters on the tray – if applicable.



7 Verify: Gently remove the FlexRx™ – 28-Cavity Weekly Standard from the tray. Your blister packs are now complete and ready for verification before dispensing to your patient.

Accessing Medication




1. Locate the appropriate medication compartment for dosage and time.



2. From the label side, push through the label of the compartment you've chosen.



3. Pinch the label and peel it back to access the medication compartment.

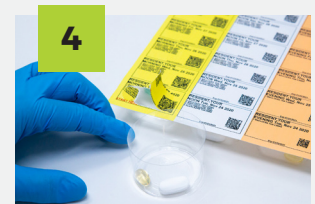

4. The medication is ready to be taken or served in a medication cup